

FrameStrip®

2-Component PCR Tube Strips

FrameStrip® 2-Component PCR Tube Strips

Strips of 8, 0.2 ml clear polypropylene PCR tubes held in a rigid polycarbonate frame.

By moulding the frame portion in a more rigid polymer, the mechanical stability of FrameStrips® is greatly improved compared with traditional single piece products, as seen in the figure 15.

For colour coding of experiments, we offer strips with clear wells in 6 different frame colours. Additionally, for optical assays such as qPCR, we supply strips with white wells in black frames. White wells increase the signal-to-noise ratio by maximising reflection of light in fluorescent based assays.

FrameStrips® are available with either domed cap strips or flat, optically clear cap strips, and are compatible with the majority of thermal cyclers. End tabs allow for easy handling and labelling of the strips and some products are also available with an off-the-shelf 2D code, offering a vast supply of unique code combinations.



Figure 15: Unlike standard tube strips, the FrameStrip® will remain straight and stable, even at elevated temperatures and when filled with liquid.

FrameStrip® 8 Well PCR Tube Strips




Features


- 2-component design in an 8 well strip format
- Compatible with the majority of thermal cyclers
- Available with either domed cap strips or with flat, optically clear cap strips
- Available with off-the-shelf 2D code

FrameStrip® 8 Well PCR Tube Strips

Standard profile, 0.2 ml PP wells, PC frame, available with either strips of domed or flat optical caps

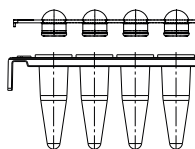

FrameStrip® 8 Well PCR Tube Strips, Plus Strips of Domed Caps

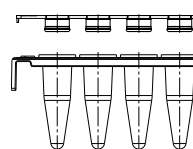

FrameStrip® 8 Well PCR Tube Strips, Plus Strips of Flat Optical Caps

Code FrameStrip® + Domed Caps	Code FrameStrip® + Flat Caps	Details		Tube and Cap Strips/Case	
4ti-0785/P	4ti-0786/P	purple frame	C	clear wells	120
4ti-0785/B	4ti-0786/B	blue frame	C	clear wells	120
4ti-0785/G	4ti-0786/G	green frame	C	clear wells	120
4ti-0785/R	4ti-0786/R	red frame	C	clear wells	120
4ti-0785/X	4ti-0786/X	black frame	C	clear wells	120
please enquire	4ti-0786/XW	black frame	W	white wells	120
4ti-0785/M	4ti-0786/M	mixed frame colours	C	clear wells	120
Code FrameStrip® 2D coded + Domed Caps	Code FrameStrip® 2D coded + Flat Caps	Details		Tube and Cap Strips/Case	
4ti-0785/X/2D	4ti-0786/X/2D	black frame	C	clear wells	120
4ti-0785/W/2D	4ti-0786/W/2D	white frame	C	clear wells	120

FrameStrips® feature a rigid polycarbonate frame for highest mechanical stability